# PROFI CARE®

## Bedienungsanleitung / Garantie

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## Instruction Manual

Thank you for choosing our product. We hope you will enjoy using the appliance.

### Symbols in these Instruction Manual

Important information for your safety is specially marked. It is essential to comply with these instructions in order to avoid accidents and prevent damage to the appliance:

#### **WARNING:**

This warns you of dangers to your health and indicates possible injury risks.

#### **CAUTION:**

This refers to possible hazards to the appliance or other objects.

**NOTE:** This highlights tips and information.

## Contents

Overview of the Components .....	3
General Notes .....	81
Special Instructions for this Appliance .....	82
Special Safety Instructions for this Appliance .....	82
Intended Use .....	85
Possible Side Effects .....	86
Hair Colour and Skin Type Table .....	88
Long-term Hair Removal .....	89
When You Should Not Use the Appliance .....	90
Unpacking the Appliance .....	90
Overview of the Components / Scope of Delivery .....	90
Application Instructions for Safe Use .....	90
What You May Expect when Using the Appliance .....	91
Before Use .....	92
Use of the Appliance .....	92
After Use .....	93

Cleaning .....	94
Troubleshooting .....	95
Technical Data .....	96
Disposal .....	97
Product Specifications in Accordance with Regulation (EU) 2019/1782 Power Supplies .....	97

## General Notes

Read the operating instructions carefully before putting the appliance into operation and keep the instructions including the warranty, the receipt and, if possible, the box with the internal packing. If you give this appliance to other people, please also pass on the operating instructions.

- The appliance is designed exclusively for private use and for the envisaged purpose. This appliance is not fit for commercial use.
- The power adapter is only suitable for use in dry rooms.
- Check if the output current, voltage and the polarity of the power adapter match the information on the connected appliance!
- Do not touch the connected power adapter with wet hands!
- Ensure that the cable of the power adapter does not become kinked, jammed, driven over or come into contact with heat sources.
- Do not use a damaged power adapter.
- Do not repair the appliance yourself, but visit an authorized specialist.
- The appliance and the mains lead have to be checked regularly for signs of damage. If damage is found the appliance must not be used.
- Use only original spare parts.

- In order to ensure your children's safety, please keep all packaging (plastic bags, boxes, polystyrene etc.) out of their reach.

**WARNING:**

Do not allow small children to play with the foil.  
**There is a danger of suffocation!**

### Special Instructions for this Appliance

This appliance is a powerful electrical appliance and, as such, should be used with special attention to the safety information provided.

Please read all warnings and safety instructions carefully before use. Precisely follow the instructions during operation as well!

The appliance is not suitable for everybody. Please follow the sections below prior to use.

A hair removal by laser or hair removal systems with an intense light pulse can cause increased hair growth in some people.

Based upon currently available data it turns out that women from the Mediterranean area, the Middle East and Southern Asia represent the largest risk group associated with side effects.

The appliance works best with darker hair types or hair that contains more melanin. Black and dark brown respond best. And although brown and light brown hair will also respond, it typically requires more sessions. Natural red hair may not respond at all. Natural white, grey or blond hair usually does not respond to applications with this appliance.

The hair growth **does not stop immediately** after the first use. You will get a permanent and satisfactory result by following the instructions in this manual.

## Special Safety Instructions for this Appliance

### Symbols on the Product

On the products you will find symbols that indicate warnings or provide information:

**Read the instructions!**

### **WARNING: Optical Radiation!**

The appliance emits strong and invisible optical radiation that may cause eye injury.

## WARNING:

Avoid irradiation of the eyes.

This appliance is **not** to be used on children.

## CAUTION:

Keep the appliance dry.

## WARNING:

Failure to follow the instructions may lead to possible eye injury, loss of vision or skin damage.

## WARNING OF ELECTRICAL SHOCK!

- Do not use this appliance near baths, wash basins or other vessels containing water.
  - If the connection cable of this appliance has been damaged, it must be replaced with a special connection cable available from the manufacturer or its customer service.
  - Keep the power adapter dry.
- Do not use any power adapter other than the one supplied.

- This appliance shall not be used by children.
- Children shall not play with the appliance.
- Appliances can be operated by persons with reduced physical, sensory, or mental capabilities or lack of experience and/or knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understood the hazards involved.
- Cleaning and user maintenance shall not be made by children without supervision.
- Do not use the appliance on areas where the skin integrity is affected. For example on skin areas with burns, injuries, blisters, scars, cuts, open wounds, active skin disorders, recent suntan or sunburn and/or infections.
- This appliance is **not** a medical appliance. Medical advice should be sought when experiencing rashes, itchy skin, cutaneous fungus and/or skin infections, skin bumps or skin spots.
- Prolonged or repeated application on one location may lead to potential damage to the skin.
- Avoid excessive passes, frequent pulses and/or increased use of the appliance.

- The appliance and light aperture optics must be regularly inspected for damage. Do **not** use a damaged appliance.
- Protect the eyes from irradiation! The appliance must not be used on the eyelids.
- A misuse of the appliance may lead to eye injury.
- Do not tamper with the built-in safety mechanism.
- Perform a skin compatibility test prior to full application. See chapter “Before Use” ⇨ “Performing a Skin Compatibility Test”.
- Please refer to the table on the different colour pigmentations in chapter “Hair Colour and Skin Type Table”. If you are still uncertain whether the use of the appliance is safe, consult a doctor or dermatologist.

## Intended Use

The appliance is used to remove undesirable body hair. The appliance is **not** suitable for individuals under the age of 18.

It is also used for reducing the regrowth of the body hair.

The appliance is designed for the following applications areas:  
Legs, arms, bikini zone, armpits, chest and back.

It is **not** suitable for the region of the nipples, anus or genitals (male or female).

## Possible Side Effects

If the instructions in this manual are followed, side effects and complications in relation to the use of the appliance will be rare. Nevertheless, every cosmetic method, including those that are intended for domestic use, is associated with certain risks. It is therefore important that you understand and accept the risks and complications associated with intense light pulse hair removal systems for domestic use.

Complication	Likelihood
Effect of multiple pulses at the same location	Minor

Side Effect	Likelihood
Slight irritation of the skin	Minor
Erythema	Minor
Increased skin sensitivity	Minor
Skin wounds and burns	Rarely
Scar formation	Rarely

Side Effect	Likelihood
Pigment changes	Rarely
Excessive redness and swelling	Rarely
Infection	Negligible
Haematoma	Negligible

### Slight Irritation of the Skin

Although a hair removal with intense light pulses for domestic use is generally well tolerated, most users experience a minor discomfort during application. The irritation of the skin is usually described as a slightly stinging sensation to the treated skin surfaces. The stinging sensation usually persists as long as the flash of light is active or even a few minutes after that. Everything that goes beyond this minor discomfort is not normal and means you should either stop using the appliance since you are intolerant to this hair removal or the selected power setting is too high.

### Erythema

Your skin may become red immediately after using the appliance or within 24 hours after the application. An erythema usually fades away within 24 hours. Consult a doctor if the erythema does not disappear within 2 to 3 days.

### Increased Skin Sensitivity

The skin of the treated area is more sensitive than usual so that it may result in dehydration or peeling off skin.

### Skin Wounds and Burns

In very rare cases, burns or skin wounds may occur after using the appliance. The burn or skin wound may take a few weeks to heal. In extremely rare cases, it will leave a scar.

### Scar Formation

Despite its rarity, however, a permanent scar formation may occur. In the event of a scar formation, it usually appears in form of a flat and white lesion (Hypotrophy) on the skin. However, it can also be large and red (Hypertrophy) or large and beyond the edges of said injury (Keloid). Subsequent aesthetic treatments may be necessary to improve the appearance of the scar.

### Pigment Changes

The appliance targets the hair shaft, in particular the pigmented cells in the hair follicles, and the hair follicles themselves. There is, however, a danger of hyperpigmentation (reinforced pigment or brown colouration) or hypopigmentation (fading) of the surrounding skin. The risk of changes in the skin pigmentation is higher for people with a darker complexion. Normally, colourations or pigment changes of the skin are temporary and a permanent hyperpigmentation or hypopigmentation occurs rarely.

### Excessive Redness and Swelling

On rare occasions, the treated skin may become very red and start to swell. This occurs more often at sensitive body sites. The erythema and swelling should subside within 2 to 7 days. Frequent cooling with ice may help to reduce the swelling. Gentle cleaning of the skin is permitted, but you should avoid exposing yourself to the sun.

### Infection

An infection of the skin is extremely rare. However, there is a potential risk subsequently to a burn or skin wound caused by the appliance.

### Haematoma

In very rare cases, blue marks (Haematoma) which can last 5 to 10 days may be formed after using the appliance. A russet coloration of the skin (Hyperpigmentation) may occur when the haematoma fades away.

## Hair Colour and Skin Type Table

Do not use the appliance on naturally dark skin or tanned skin. The application on dark skin may result in burns, blisters and changes in skin colour.

Hair Colour	Celtic Skin Type	Nordic Skin Type	Mixed Type	Mediterranean Skin Type	Dark Skin Types	Black Skin Types
	I	II	III	IV	V	VI
Black	4/5*	4/5*	3/4*	1/2/3*	X	X
Dark brown	4/5*	4/5*	3/4*	1/2/3*	X	X
Brown	4/5*	4/5*	3/4*	1/2/3*	X	X
Dark blond	4/5*	4/5*	3/4*	1/2/3*	X	X
Blond	X	X	X	X	X	X
Red	X	X	X	X	X	X
Grey & white	X	X	X	X	X	X

\* = applicable light intensity level X = not suitable

## Long-term Hair Removal

A typical hair growth cycle can last between 18 and 24 months. In order to achieve a long-term hair removal, multiple sessions with this appliance may be necessary during this period. Depending on hair colour, the body site and how the appliance is being used, the effectiveness of the hair removal will vary from person to person.

## Typical Application Plan during a Hair Growth Cycle:

Number of Sessions	Time Interval in Weeks
1-4	2
5-7	4
8+	As needed, until the desired results have been achieved

## The Expected Results:

Fitzpatrick Skin Type	Hair Colour	Anatomical Body Region	Number of Sessions	Effectiveness*
I-IV	Light brown to brown	Upper or lower legs	Around 10 to 12 with 4 weeks between the sessions	60 % Decrease of the hair growth
		Arms		50 % Decrease of the hair growth
		Bikini zone		60 % Decrease of the hair growth
		Armpit		55 % Decrease of the hair growth
I-IV	Brown to black	Upper or lower legs	Around 8 to 10 with 4 weeks between the sessions	60 % Decrease of the hair growth
		Arms		50 % Decrease of the hair growth
		Bikini zone		60 % Decrease of the hair growth
		Armpit		55 % Decrease of the hair growth

\* (% of hair removal/renewed hair growth)

- The response is individual, depends on the hair type and biological factors, and may affect renewed hair growth. Some people show a result sooner or only at a later time.
- After one to two sessions a permanent decrease of hair growth cannot yet be expected. Furthermore, the inactivity of the hair follicles varies depending on the body site.

## When You Should Not Use the Appliance

Do not use the appliance if any of the following statements applies to you:

- If you are pregnant or breastfeeding.
- If you have been exposed to intensive solar radiation for the last 4 weeks or if you have been to a solarium.
- In areas with tattoos, permanent make-up or over piercings and other metallic items.
- In areas with large freckles, birthmarks, moles or warts.
- If you are suffering or have suffered from keloid scars.
- If you are sensitive to light or on medication that makes the skin more sensitive to light.
- If you are suffering from abnormal skin conditions because of diabetes, for instance, or if you are suffering from another metabolic disorder.
- If you recently had a treatment with, for example, alpha hydroxy acids (AHAs), beta hydroxy acids (BHAs), Retin-A®, Accutane®, topical retinoids and / or azelaic acid.
- If you have received a steroid treatment in the last 3 months.
- If you have a history of herpes in the treatment area.
- If you suffer from epilepsy.
- If you have or have been fallen ill with skin cancer.
- If you have an active implant such as a pacemaker, an appliance for incontinence, an insulin pump, etc.
- If you have undergone a radiation therapy or chemotherapy in the last 3 months.
- If you are suffering from a disorder of the blood coagulation.
- If you have other diseases that, according to your doctor, would make the treatment unsafe for you.

If you are still uncertain whether the use of the appliance is safe, consult a doctor or dermatologist.

## Unpacking the Appliance

1. Take the appliance out of its packaging.
2. Remove all packaging materials, such as foils, filler material, cable brackets and cardboard packaging.
3. Check the contents for completeness.
4. In case of incomplete delivery or damage, do not operate the appliance. Return it to your dealer immediately.

## Overview of the Components / Scope of Delivery

- 1 Air outlet openings
- 2 Connection socket for power adapter
- 3 Indicator lamp START
- 4 ⏻ button (on / off / select light intensity)
- 5 Indicator lamps and light intensity level display
- 6 Release button for the light flash
- 7 Light flash outlet window
- 8 Skin contact surface
- 9 Attachment
- 10 Power adapter

### No illustration

Safety glasses

## Application Instructions for Safe Use

- **Carefully Select the Intensity Level of the Light Flashes**

The intensity level refers to the intensity of the emitted light flashes during use, from the lowest to the highest level. The indicator lamps next to the level display (5) show the intensity level that has been set. Higher intensity levels improve the results, but also the resulting risks and side effects. (See page 82 "Possible Side Effects")

- **Always Start the Application on the Lowest Intensity Level**

Only increase the intensity level for the next application by one step if you experience little or no discomfort during or after use of the appliance while having set the lowest intensity level. However, also note the indicated intensity level in the table on hair colours and skin types.

- **Avoiding Side Effects**

- Do not use the appliance multiple times at the same skin area during one session!
- Avoid overlapping light flashes!
- Immediately stop the usage when your skin shows blisters or burns!

- **Avoid Complications After Use**

- Do not expose the skin to the sun after the application. This also applies to tanning with artificial light (Solarium).

The exposure to the sun involves a continuously unprotected solar radiation for 15 minutes. And a constantly unprotected solar radiation of diffuse sunlight for more than one hour. In order to protect the treated skin against sunlight, apply a sunscreen with a sun protection factor that is equal or higher than SPF 30 for 2 weeks after each use of the appliance.

- Do not use the appliance within 24 hours before or after visiting a sauna or swimming pool. The chlorine could lead to skin irritation.
- Cover birthmarks and tattoos before starting the application.
- Cover dark brown or black spots such as freckles, moles or warts before starting the application.
- Never look directly into the light flash! **Always** wear the supplied safety glasses when applying.
- Do not use the appliance in the region of the nipples, anus or genitals (male or female).
- Do not use the appliance at body sites where you might want to be hairy later.

- Do not use the appliance for any other than the intended purpose.
- Always ensure that the skin contact surface is in full contact with the skin.
- Remove the appliance from the skin if either the skin contact surface or the skin feels hot.

## What You May Expect when Using the Appliance

Perhaps for many users it is the first experience with a light-based hair removal appliance that has been developed for domestic use. The appliance is easy to use and hair removal sessions are hardly time-consuming. It is normal that you perceive the following during a session with this appliance:

- **Light Flash**

A light flash may be visible when using the appliance. Therefore, **always** wear the supplied safety glasses when applying.

- **Noise of the Fan**

The built-in cooling fan of the appliance makes a sound that is similar to a hair dryer. This is normal.

- **Popping Sound with Every Light Flash**

It is normal that you hear a popping sound at the same time as the light flash is triggered.

- **Slight Imprint of the Light Flash Outlet Window**

This is necessary and useful, and a part of the appliance's safety feature to ensure the placement of the adjacent light flashes.

- **A Sense of Warmth and Tingling**

During each light flash it is normal that you experience a slight feeling of warmth and tingling. Remember that it is important to always use a low power level for the first hair removal sessions. There may be times when you feel certain warmth up to one hour after the application.

- **Slightly Red or Pink Skin Colour**

During or directly after a session, it is not uncommon to note a very slightly pinkish colouration of the skin. This pink colour is usually visible around the region of the hair itself. However, if your skin is completely red or blistering and burns occur, stop using the appliance immediately.

## Before Use

### WARNING:

**Never** use flammable liquids such as alcohol (including fragrances, disinfectants or other alcoholic substances) or acetone to clean the skin before using the appliance.

- Shave the skin at the location you want to use the appliance. The hair must not be removed with wax or ripped out in any kind. The shave is an important prerequisite for the desired result.
- For best results, shave the region 3 days prior to application or when the hair has a length of 1 to 2 mm.
- The skin should be clean, dry and free of any antiperspirant or deodorant.

### Performing a Skin Compatibility Test

The skin compatibility test is designed to determine how your skin responds to the treatment with the appliance. Furthermore it aids in adjusting the optimal power level preference for each treatment area.

Perform a skin compatibility test on each skin area you would like to treat **48 hours** before the full treatment.

### WARNING:

**Always** wear the safety glasses when performing a skin compatibility test!

1. Proceed as described in chapter “Use of the Appliance” point 1 to 4.

2. Start with the lowest intensity level.
3. Place the skin support surface onto the desired test location on the skin.
4. Proceed as described in chapter “Use of the Appliance” point 6 to 8.
5. Continue the test with the next higher intensity level on another test location on the skin.  
However, also note the indicated intensity level in the table on hair colours and skin types.

- Repeat the skin compatibility test for each body region in which you want to use the appliance.
- If you don't experience any discomfort within 48 hours, you may start the application with the last selected intensity level.

## Use of the Appliance

### WARNING:

- Perform a skin compatibility test before using the appliance.
- **Always** wear the safety glasses when applying!

### NOTE:

The appliance is equipped with a built-in safety mechanism. This prevents, if the skin contact surface is not fully placed on the skin, the generation of a light flash when pressing the release button.

1. Make sure the grid voltage corresponds to the figures on the rating plate.
2. Plug the power adapter unit into a properly installed protective contact socket.
3. Connect the appliance with the power adapter unit by plugging the jack of the power adapter unit into the jack (2). All indicator lamps for the level displays light up shortly. The indicator lamp START lights red.

 **CAUTION:**

- This appliance should only be used with the power adapter unit (included). Do not use it for other appliances.
  - Only use a 12 V power adapter unit (⊖—C—⊕). A different power adapter unit may damage the appliance.
4. Press the  button to switch the appliance on.
    - The cooling fan of the appliance will start.
    - The indicator lamp of level display 1 is lit.
    - And the indicator lamp START switches from red to green after approx. 1 second. The appliance is now ready for use.
  5. To change the intensity level, press the  button repeatedly. The corresponding indicator lamp of the level display is lit.

**Intensity level**

The intensity level determines the intensity of the emitted light flash. There are 5 levels.

1 = lowest level, 5 = highest level

The indicator lamp next to the level display shows the intensity level that has been selected. With a higher intensity level the hair removal result will improve, but also the risk of possible side effects and complications will increase.

Do not increase the intensity level within a session. Symptoms may only occur after a couple of hours, even if you do not experience any of them at first.

6. Place the skin contact surface onto the skin location. Make sure that the skin location is even and flat, and that the skin contact surface makes full contact with the skin.
7. Once the skin contact surface rests completely on the skin, the indicator lamp START will start flashing. The appliance is now ready to emit a light flash.

 **NOTE:**

The indicator lamp START only starts flashing when the skin contact surface fully rests on the skin.

8. Press the release button for the light flash. The appliance will emit a light flash and you will hear a short popping sound. This is normal for the appliance. At the same time you may experience a mild, warm and tingling feeling.
9. After the appliance has emitted a light flash it will power up for the next light flash. This may take a few seconds.
10. Place the skin contact surface onto the next skin location. Make sure to irradiate the same skin location only once per session.
11. Press the release button for the light flash.
12. Proceed as described from point 6.


 **WARNING:**

**Immediately stop the usage if your skin forms blisters or starts to burn!**

 **NOTE:**

The appliance turns off automatically after 10 minutes.

## After Use

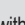1. Turn off the appliance when the session is done. To do this, press and hold the  button for about 3 seconds. Remember the last selected intensity level. The appliance generally starts with the lowest intensity level the next time you turn it on.
2. Disconnect both the power adapter from the socket and the power plug from the appliance.
3. Clean the appliance and especially the skin support surface after each session. (See chapter “Cleaning”)

### Skin Care After the Hair Removal Session

Do not expose the treated skin to the sun. Make sure to protect the treated skin against sunlight. Apply a sunscreen with a sun protection factor that is equal or higher than SPF 30 for 2 weeks after each use of the appliance.

### Side Effects and Complications

Some users may experience pigment changes due to the treatment with the appliance. These changes are, if they occur, generally mild and temporarily. If complications occur, contact your dermatologist immediately!

## Cleaning

### WARNING:

Switch off the appliance before cleaning. Disconnect the power adapter from the electrical outlet.

### CAUTION:

- Do **not** immerse the appliance or any other part of it in water!
- Do not use a wire brush or other abrasives.
- Do not use sharp or abrasive cleaning agents.

Use a dry, clean cloth and, if necessary, a cleaner that is designed specifically for electronic equipment to gently wipe the surface of the appliance and in particular the skin contact surface.

### Storage and Maintenance

- Keep the appliance in the original packaging after cleaning. Store the appliance in a cool and dry location between each session. Keep the appliance away from water.

- Regularly check the light flash outlet window for damage or scratches. Do not operate a damaged appliance!
- The attachment needs to be replaced after approx. 50000 light flashes.

You can order spare attachments from your local dealer. (Spare parts no.: 399995)

## Replacing the Spare Cartridge

You can easily replace the attachment yourself.

1. Remove it by simply pulling it forwards.
2. Slide the new attachment onto the appliance.

Ensure that the attachment has been fitted properly, as otherwise the appliance will not emit a light flash.

### Troubleshooting

Problem	Possible Cause	Solution
The appliance does not function.	The appliance has no power supply.	Ensure that the connector of the power adapter is properly plugged into the appliance.
		Ensure that the power adapter is properly plugged into the socket.
		Test the socket with another appliance.
		Check the main circuit breaker.
	The appliance is faulty.	Contact our service centre or a specialist.

Problem	Possible Cause	Solution
It will not emit a light flash when the release button is pressed.	The attachment is not properly placed on the appliance.	Ensure that the attachment has been mounted correctly.
	The appliance is not yet operational.	Makes sure that the indicator lamp START flashes green.
	Poor contact with the skin.	Ensure that the skin support contact fully rests on the skin.
	The lamp in the attachment is defective.	Replace the attachment with a new one. If the error persists, contact our customer service.
Smoke can be seen or smelled during operation.	The appliance is faulty. Immediately stop the usage!	Contact our service centre or a specialist.
The indicator lamp START and one of the level displays remain lit.	System error	Switch the appliance off. Disconnect the power adapter from the socket for around 5 seconds. Re-insert the power adapter into the socket and switch the appliance back on. If the error persists, contact our customer service.
The indicator lamp START is lit. And the other indicator lamps start flashing.	The appliance has reached the end of its life.	Dispose the appliance properly. See chapter "Disposal".

## Technical Data

Model:.....PC-IPL 3024  
 Power supply:.....DC 12V, 2,5A  
 Light surface:.....2.7 cm<sup>2</sup>  
 Technology used:.....IPL (Intense Pulsed Light)  
 Flash duration: .....500 to 800 µs  
 Intensity density: .....4J/cm<sup>2</sup>  
 Wavelength: .....475 to 1200 nm  
 Net weight: .....approx. 0.4 kg  
 Temperature:  
   Operation:.....+10 °C to +35 °C  
   Storage:.....-40 °C to +70 °C

Relative air humidity:

Operation:.....3 % to 75 % rH  
 Storage:.....10 % to 90 % rH

Atmospheric pressure:

Operation:.....700 to 1060 hPa  
 Storage:.....500 to 1060 hPa

## External Power Adapter

Input: .....AC 100 – 240V, 50/60 Hz, 0,8 A Max  
 Output:.....DC 12V/3A  
 Polarity: .....⊖—⊕  
 Protection class:.....II

The right to make technical and design modifications in the course of continuous product development remains reserved.

This appliance has been tested according to all relevant current CE guidelines, such as electromagnetic compatibility and low voltage directives, and has been constructed in accordance with the latest safety regulations.

## Disposal

### Environmental Protection

Do not dispose of the appliance at the end of its service life with normal household waste. Bring it to the official collection point for recycling. By doing this, you help protect the environment.

## Product Specifications in Accordance with Regulation (EU) 2019/1782 Power Supplies

In case the technical data deviate from the data listed in the operating instructions, only these in the table below are valid.

Information about the power supply used	Value	Unit
Clatronic International GmbH Commercial register number HRB 9158 Industriering Ost 40, 47906 Kempen Germany		
Trade mark	ProfiCare	
Model identifier	PC-IPL 3024	
Input voltage	100 - 240	V
Input alternating current frequency	50/60	Hz
Output voltage	12.0	V
Output current	3.0	A
Output power	36.0	W
Average efficiency in operation	87.6	%
Low load efficiency (10%)	86.3	%
No-load power consumption	0.06	W

The relevant load conditions are as follows:

Percentage of nameplate output current

Load condition 1	100 % $\pm$ 2 %
Load condition 2	75 % $\pm$ 2 %
Load condition 3	50 % $\pm$ 2 %
Load condition 4	25 % $\pm$ 2 %
Load condition 5	10 % $\pm$ 1 %
Load condition 6	0 % (no-load condition)

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Clatronic International GmbH Numer rejestru handlowego HRB 9158 Industriering Ost 40, 47906 Kempen Niemcy