Vögan<sup>®</sup>  
Mammoth

*Lady Two*

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# Vacuum skin cleaner with micro dermal abrasion **Vigan Mammoth Lady TWO**

Dear customer, we would like to thank you for buying Vigan Mammoth Lady TWO (further only "device"). Please read carefully this instruction manual before usage. Use this product only as described in this manual. Keep the instruction manual for later use. We recommend to keep the original packaging and receipt. If the device will be handed to another person, always attach the instruction manual.


## GENERAL INFORMATION

Vacuum microdermabrasion is one of the most popular aesthetic dermatology treatments. Thankful to this device, you can enjoy this procedure from your home. This device cleans and rejuvenates your skin and reduces facial wrinkles. This method effectively helps to remove the black dots (comedones) from your nose and eliminate the skin pigmentation.

The special head with fine microscopic crystals (is included in package) removes the micro layers of dead skin cells from the skin's epidermis. This treatment enhances cell renewal and regeneration of the skin surface. The treatment significantly improves the skin's absorption. The continuously created vacuum improves blood circulation and combination of all effects gradually leads to increased production of collagen and elastin fibers. The result is the strengthened layers of your skin which improves overall skin health.

Microdermabrasion does not cause any major contradictions, its harmless and suitable for almost every skin type. This treatments features non-invasiveness and zero risk.

Vacuum lifting and peeling performed by Lady TWO works on a vacuum principle which gently sucks in the skin and due to under-pressure the skin is cleansed and stretched. You can choose from 5 different attachments that are included in the package.

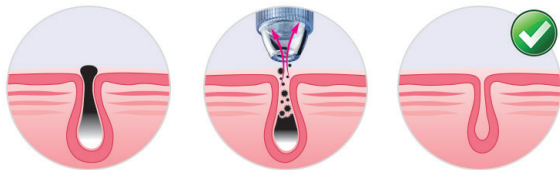Lady TWO, USB cable, 5 attachments, 50pcs of spare micro-filters, 2pcs of silicone sealing rings

**Vigan**  
Mammoli

str 31

## 5 TYPES OF EXTENSIONS FOR DIFFERENT USES

### 1. Small round extension

- Suitable for weak and sensitive skin
- Suitable for low-suction mode

### 2. Large round extension

- Suitable for effective blackhead removal
- Suitable for skin cleansing
- Cleans skin with slow movements

### 3. Oval extension

- Suitable for badly reachable face creases

### 4. Micro-crystal peeling extension

- Removes dead skin
- Enhances cleansing
- Stimulates rejuvenation

### 5. Blue light extension


- This extension has antiseptic and anti- inflammation effect
- Calms down the skin
- This extension is used at the end of the treatment

### Micro filter interchange:

**The key for the right treatment is to open the pores**  
the pores can be opened by hot shower, spa or steam
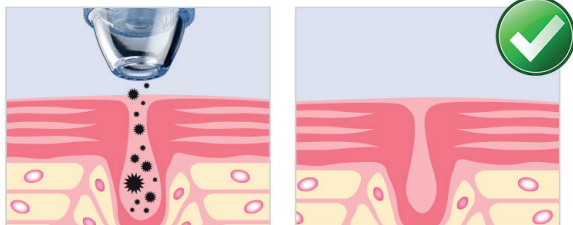
## MAIN CHARACTERISTIC OF VIGAN LADY TWO

1. Device is designed for easy and effective skin care
2. 3 levels of suction intensity
3. Blue light therapy – calms the skin and has antiseptic and anti-inflammation effect
4. Non-invasive, painless treatment for almost every type of skin
5. Vacuum lifting and peeling from the comfort of your home
6. Improves blood circulation and has anti-aging effect
7. Removes dead cells and supports its regeneration
8. Removes the impurities in the pores
9. Supports the microcirculation and strengthen the layers of your skin
10. Device stimulates the collagen and elastin production
11. Ergonomic shape for easy handling and manipulation
12. Easy to carry
13. Low weight and space saving design
14. USB charging
15. Up to 60 minutes of operation time

## **Safety instructions:**

- Search the device and accessories for any damages. If you find any, don't use the device and contact the seller or professional after sale service.
- Never immerse the device into the water! Do not use the device with wet hands.
- Do not use other accessories than the original delivered by the manufacturer.
- Do not leave the device unattended while turned on.
- Dispose the packaging materials in appropriate way. Beware! The plastic bags can be dangerous. To avoid risk of suffocation, keep all the plastic waste away from children and babies.
- Use and store the device in dry environment!
- Do not charge the device close to bath, spa, pool etc.
- Never use the device if it fell into the water.
- Protect the device from shocks, dust, chemicals and sources of heat.
- Be careful not to drop the device.
- Make sure that the charging cable will never come into contact with water/ steam or other fluids.
- Device is not suitable for users of age under 18 years.
- Device is for external use only.
- This device is designed for human face treatment (except eyes, mouth, ears).
- Never use on animals.
- **This device can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they are properly instructed and supervised about safe and proper use of this equipment. Cleaning and user maintenance by the children should not be done without supervision. Children must not play with the device.**
- If the charging cable is broken, it has to be replaced with the same type of cable, otherwise the device can be damaged and become dangerous.
- Don't use the attachments if they are broken, worn or dirty.
- This device includes magnet. Keep the device in distance from pacemakers. If persons using pacemakers (or other implanted electronic device) are interested in using this device, they have to first discuss with the doctor!
- Persons using pacemakers must stay away from the device.
- Do not disassemble the device! Do not change the batteries!
- If there is any issue with the device, always contact the seller or professional service.
- Do not use the device for purposes other than those specified in this instruction manual. Manufacturer can not hold responsibility for damages caused by improper use.

## **Principles to avoid skin irritation:**

- Do not use Lady TWO more than twice a week.
- Do not combine the usage of peeling attachment with other attachments.
- While using, always keep the device in motion – never stop on one spot! The vacuum could cause irritation/ bruise or other injury!
- Do not combine the treatment of the device with body lotion/cream or peeling cream. This could lead to skin injury or device damage.

## **Do not use:**

- If you have cracks on your skin
- If you have reddish, irritated, damaged, bleeding skin
- If you have any skin disease/ burns
- On eyes, lips and ears
- If you suffer from skin cancer
- If you use blood thinners
- Directly after sun tanning/ If your skin is sunburnt
- If you have hypersensitive skin
- On birth marks/ tattoos/ blisters/ inflamed areas/ severe acne
- At least two weeks if you underwent cosmetic surgery
- If you are using skin care medicine/ treatments
- If you suffer from multiple sclerosis, diabetes or any autoimmune diseases

- **IF YOU HAVE ANY SKIN DISEASES, CONSULT YOUR DOCTOR BEFORE USAGE OF THE DEVICE.**
- **IF YOU DON'T FEEL COMFORTABLE WHILE USING THIS PRODUCTS/ HAVE SKIN PROBLEMS, CONSULT YOUR DOCTOR**
- **USAGE OF THIS DEVICE DOESN'T SUBSTITUTE THE MEDICAL TREATMENT - CONSULT YOUR DOCTOR BEFORE THE USAGE**

## CHARGING

**We highly recommend to fully charge the device before every use.**

- If you don't use the device for more than three months, charge it. It will extend the battery life.
- The device can be used while charging.
- Charging takes approximately 3 hours. The battery life is around 60 minutes (depends on intensity level and type of attachment used).
- Before you plug in the cable, make sure that the device is turned off.
- First plug the cable into the device and after that, connect the other end to USB charging port.
- All of the blue diodes flash during the charging.
- When the device is fully charged, the blue diodes will light.
- When the battery is fully discharged, the device will automatically turn off and the diodes will switch off.

## BEFORE USAGE

Always wash your face thoroughly and dry before the usage. Make sure that you removed the make-up, cream etc. We recommend steaming your face before the treatment- it will open your pores and improve the effect.

## AREAS OF FACE TREATMENT

**Nose:** Move the device slowly from the bridge of your nose down to its tip and sides.

**Temples:** Move the device from outward to the hair line. Be careful to keep the device at least 1 cm from your eyes!

**Forehead:** Move the device slowly from the center of the forehead towards the hair line. Make sure that you don't treat any spot more than once!

**Cheeks:** Move the device slowly from the sides of your nose towards your jaws.

### Surrounding of the mouth:

Move the device from the center under the nose/ chin to both left and right sides. Make sure that the device will never come in contact with your lips!

## FIRST USE OF THE DEVICE

We strongly recommend to conduct the device test on your forearm before the first use.

1. Attach the device head that you wish to use on your face. Make sure that there are the silicone seal rings under the attachment. You can use circular movement to install the attachment properly.
2. Turn on the device
3. Choose the desired intensity (3 levels + 4th with blue light with skin calming effect)
4. Before you start the treatment on your face, apply the device on your forearm and make couple of draws and wait approximately 30 minutes. If there is no severe reaction, you can start the treatment on your face. If there is any unexpected reaction, consult your doctor!

## FIRST FACIAL USE OF THE DEVICE

For the first use on your face, we strongly recommend to start with the lowest level of suction so your skin can get used to this treatment. If there will be no difficulties after first application, you can increase the suction power.

1. Attach the device head that you wish to use on your face. Make sure that there are the silicone seal rings under the attachment. You can use circular movement to install the attachment properly.
2. Turn on the device
3. Choose the desired intensity (3 levels + 4th with blue light with skin calming effect)
4. Apply to your skin so you can feel the suction. Make sure that the attachment is in full contact with your skin. There is no need of pushing the device against the skin too hard.
5. Try to make couple of draws on your forehead and wait 30 minutes. If there will be no reaction, you can use the device on the rest of your face.
6. Its strongly recommended to use the device twice a week with at least 48 hours period between the treatments.

## AFTER THE USAGE

After the usage of the device apply face mask or hydration cream.

**At least three days after the treatment obey following rules:**

- Never use any mechanical or chemical peeling!
- Do not attend any other facial treatments!
- Protect your skin from sun rays!

# CLEANING AND MAINTENANCE

## **CAUTION:** Never immerse the device in the water!

Never use rough wipes, abrasive scrubbers or aggressive chemical solution to clean the appliance!

### Cleaning:

1. Make sure that the device is turned off.
2. Carefully remove the attachment.
3. Clean the attachment with hot water. You can use the dishwashing detergent to remove the grease from the surface of the attachment.
4. Immediately after cleaning dry the attachment with the cloth and put it back into the case.
5. Change the attachment filter if necessary.
6. Clean the device with wet cloth, dry it and store it into the case.

# PROBLEM SOLVING

Please find below most common questions and answers that can occur during usage.

## Most common problems and their solving

### **Is this device suitable for my skin?**

Lady TWO is suitable for almost all kinds of skin. Please find the exceptions of the usage in the general information.

### **How long does the treatment take?**

The length of the treatment may vary but the complete face treatment should take around 20 minutes.

### **Is this technology safe?**

This device uses microdermabrasion principle, which removes the unnecessary skin particles. It is non-invasive and safe. To make it as safe as possible, read the safety instruction in this manual.

### **The device doesn't work?**

Try to charge the device as described in this manual, eventually try to change the USB cable/ source of power. If the problem persists, contact the after sale service or the importer directly.

### **Is the Lady TWO treatment painful?**

The treatment should be painless. Application on highly sensitive skin can cause red spots which will disappear. If the treatment is painful, lower the level of the suction/ try to change the attachment. If the problems persist, consult your doctor. To avoid any kind of bruises/ swellings, make sure that you treat every spot on your face only once. All these instructions are described in the instruction manual and should be read carefully before the usage.

### **Can I sunbath after the treatment?**

You should avoid the sun beams for at least 3 days.

### **The device movement is not steady/ there is weak or no suction/ the results of the treatment are not satisfying?**

Check the microfilter/ attachment if its not dirty or clogged. If the filter is dirty, change it. if the attachment is dirty/ clogged, clean it. Always use the procedures as described in „Cleaning and maintenance“ paragraph. You can also check the suction performance by laying your finger on the device without the attachment. If you wont feel any suction power, try to charge the device properly. If it won't help, contact the professional service.

### **My skin and its reaction to the treatment?**

Right after the treatment your skin may become red. During the treatment, you can feel the warmth caused by the improved blood circulation and tingling. These effects will disappear in couple minutes. Your skin can also appear swollen, red, dry and full of skin flakes. These effects may be caused by inappropriate pressure during the treatment or intensive care on one spot. To reduce these effects, don't push the device against your skin too hard and make sure that you do not stay on one spot for too long. The device was tested and it is safe to use. Always make sure to use it as described in this manual.

### **It's difficult to maintain steady suction.**

Try to tilt the device while usage and find the ideal angle of the attachment so the contact between the device and the skin will not be interrupted. While treating the forehead area, move the device slowly.

### **Can i perform any other face treatments after the device usage?**

Do not undergo any mechanical or chemical peelings, electrolysis or face epilation in following 3 days after the LadyTWO treatment.

## TECHNICAL SPECIFICATION

Battery type/capacity: Lithium battery 600 mAh

Charging voltage: 5,0 V

Operating voltage: 3,7 V

Power output: 0,8 W

Charging time: approximately 3 hours

Operation time: up to 60 minutes

Gross weight: 381g

Net weight: 135 g

Gross dimensions: 177 x 55 x 215 mm

Net dimensions: 148 x 55 x 50 mm

## PACKAGE INCLUDES

1x Lady TWO device

1x Large circular attachment

1x Peeling attachment with micro crystals

1x Oval attachment

1x Blue light therapy attachment

1x USB cable

2x Silicone seal ring

50x Spare micro filter

1x Instruction manual (CZ/ SK/ UK/ DE/ HU)

# DISPOSAL OF OLD ELECTRICAL AND ELECTRONIC EQUIPMENT

This symbol on the product, its accessories, packaging or accompanying documents indicates that the product should not be treated as household waste. Please take this product to the respective collection point where will be carried out specialist recycling of electronic equipment. In the European Union and other European countries there are collection points for used electrical and electronic equipment. By ensuring proper product disposal, you prevent potential negative impact for the environment and human health, which could result from improper handling of this product. Recycling of these materials also helps to preserve natural resources. For this reason, please do not throw the old electric and electronic equipment together with household waste. For more information please contact your local authority, any service company liquidating household waste or the store where you purchased the product. Details also refer to the relevant laws about waste, in the Czech Republic no. 185/2001 Coll. as amended. Furthermore, the websites [www.elektrowin.cz](http://www.elektrowin.cz), [www.remasytem.cz](http://www.remasytem.cz) and [www.asekol.cz](http://www.asekol.cz) (disposal of electrical waste). The improper disposal of this waste can be in accordance with applicable laws, regulations and fines.

## WARRANTY

We offer warranty for 24 months starting from the day of purchase by you. This warranty covers defect caused by defective material or error in workmanship. It does not cover defects caused by mechanical damage, natural disaster or other inevitable events, connecting to the wrong power grid voltage, modifications or service work by non-qualified or unauthorized personal, misuse or use in contradiction to this manual.

This device is equipped with long-life Li-ion battery. The 24 months warranty doesn't include the decreasing battery capacity caused either by the use or age. The capacity decreases depending on how often and how do you use/charge the battery.

All documents about your purchase and service repairs on any product keep for at least three years to ensure quality customer service.

Please remember, that this product is only for home use, not professional use.

If you find out, that this product does not suit you, it is not a reason for a complaint.

In the case of complaint please submit receipt.

In the case of malfunction of your product, please contact the dealer, you bought it from, or service at the following address:

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